



**Budget friendly** NEW **Opti-Seal™ Tear Off Optical disposable for closure of (q)PCR Strip tubes & (q)PCR plates**

For Diagnostic & Clinical procedures

- Designed to perform diagnostic qPCR and PCR reactions
- Can also be used in Real Time PCR applications
- Fits All types of PCR and qPCR plates
- Fits all types of 0.1 ml and 0.2 ml PCR and qPCR 8-tube strips

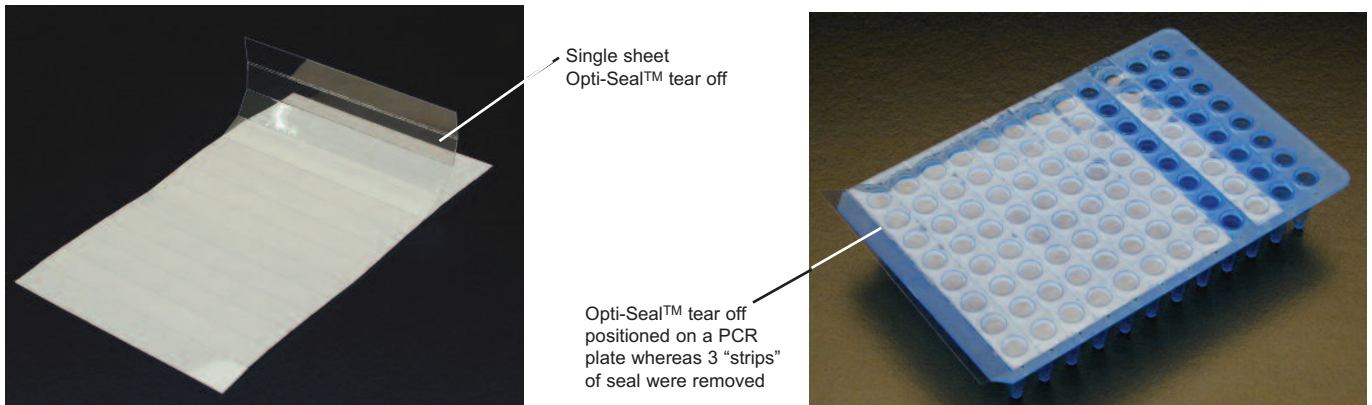

FEATURES	ADVANTAGES	BENEFITS
Tear Off capability	Individual closure of tube strips	Minimize cross contamination
		Economic usage
	Closure of plates or sections of a plate	Economic usage
	Opening of individual rows of a plate	Minimize cross contamination
Optical clear	Higher S/N ratios qPCR applications	Low Copy detection
Uniform & thick "Elastic Adhesive"	Easy Opening and Closing	"finger" and RSI friendly
	Minimize evaporation	Consistent reaction conditions
	Seal mainly based on Adhesive properties	Excellent seal for vessels with flat chimney ledge
	Seal mainly based on "rubber elastic" properties	Superior seal for vessels with non flat (angled)chimney ledge

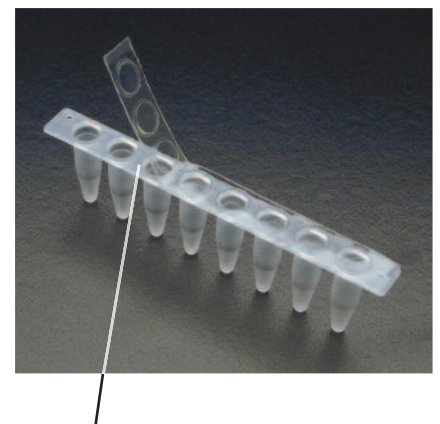

Opti-Seal™ tear off single strip

**€ 115 for 100 sheets**

**EU Opti-Seal™ Tear Off Optical disposable adhesive**

Opti-Seal™ provides the best sealing option for (p)PCR plates and (q)PCR 8-strip tubes- The EU Opti-Seal™ is non pierceable and can be easily removed after the (q)PCR reaction is performed. Pressure applied by the heated lid of the thermal cycler keeps the seal well closed during thermal cycling. Opti-Seal™ generates superior results and is designed and tested to be used in Real-Time PCR applications.


Order no	Description	Package Size	Price €
157400	EU Opti-Seal™ Tear Off Optical Adhesive	100 sheets	115

**Related products**

Order no	Description	Package Size	Price €
157300	Opti-Seal™, optical, disposable, adhesive	100 sheets	91.62
157200	EU Alu-Seal disposable, adhesive	100 sheets	91.22
B57601	EU Optical Wide Area Indented 96 cap-plate, cuttable, natural	25 plates	50.50
Miscellaneous	8, 12 and 96 cap strips	see www.bioplastics.com	





Order number	Name	Material	Adhesive	Tear Off multiple rows	Long Term Cold Storage	Pierceable	Disposable	Suitable for 384 wells	Optical for qPCR	Easy Removal
157300	EU Opti-Seal™	Polyester	Y	N	Y	N	Y	Y	Y	Y
157400	EU Opti-Seal™ Tear Off	Polyester	Y	Y	Y	N	Y	N	Y	Y
	Y=YES / N=NO									

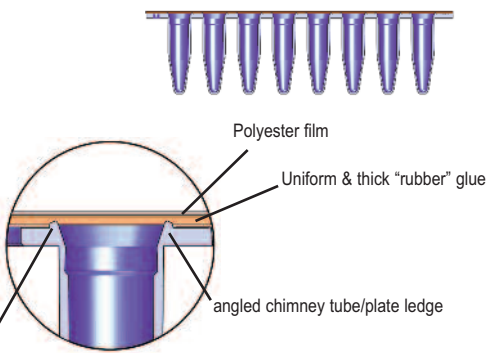


**"rubber" seal**

Opti-Seal™ seal mainly based on "rubber" seal capacity, provides excellent seal for vessels with non flat (angled)chimney ledge

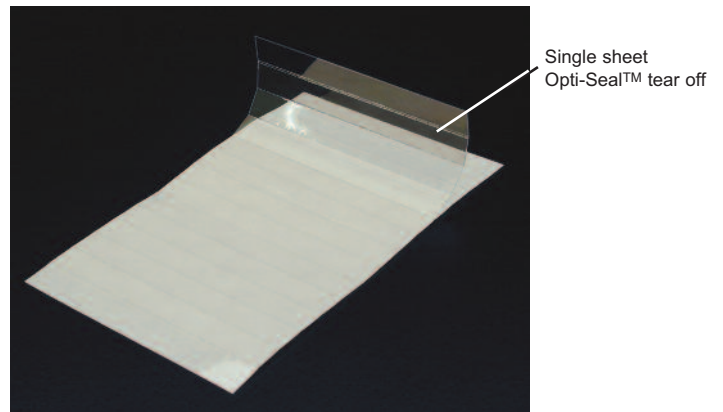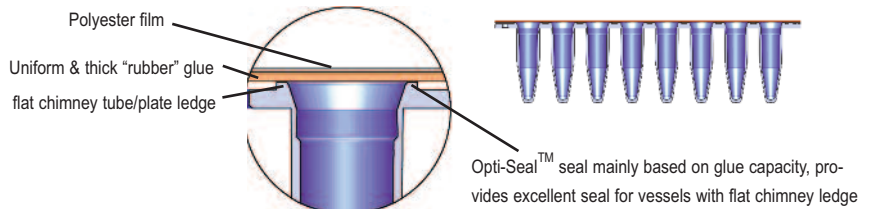

Single sheet Opti-Seal™ tear off

Opti-Seal™ tear off single strip

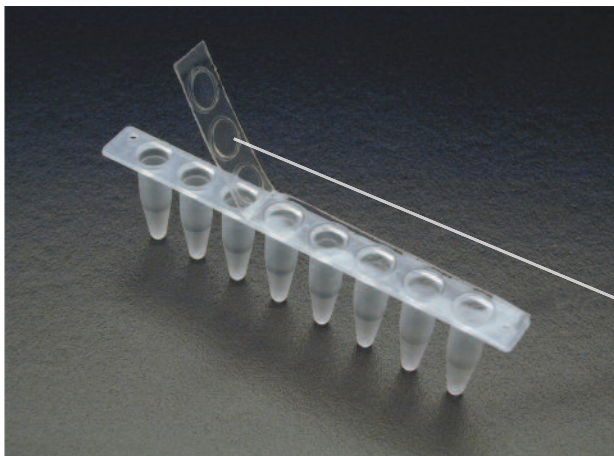

**Trademarks**

(registered) Trademarks and Owner:

Opti-Seal, Opti-Seal tear off, BIOplastics (logo), EU, BIOplastics owner BIOplastics BV

Products are designed and manufactured in the Netherlands. SAMPLES on request and/or purchase with money back guarantee\* (\* if proven to be not working) Contact your local representative. All rights reserved. 1196VKC1 ©BIOplastics BV

**Headquarters:**

(R&D and Manufacturing)  
 BIOplastics BV  
 Röscherweg 61  
 Entrance D  
 6374 XW Landgraaf  
 The Netherlands

Phone: +31 (0)45 533 87 50  
 Fax: +31 (0)45 533 87 96  
 Website: www.bioplastics.com  
 E-mail: info@bioplastics.com

**Subsidiary:**

BPCTi Inc.  
 2933 South Miami Blvd.  
 Suite 121  
 Durham, NC 27703, USA

Phone: 919 806 8811  
 Fax: 919 806 2014  
 Website: www.bpcti.com  
 E-mail: info@bpcti.com

**Your local distributor:**

